



ANALYTICAL METHOD DEVELOPMENT AND VALIDATION OF ANTIDIABETIC DRUGS (METFORMIN AND LINAGLIPTIN) IN TABLET DOSAGE FORM BY USING RP- HPLC METHOD

**Anjali Chouhan, Neelam Rathore, Anubhav Shrivastava, Mahendra Singh Tomar,
Shubham Singh Panwar**

Assistant Professor

Dr. A. p. j. Abdul Kalam University, Indore

ABSTRACT

A new, simple, precise, accurate and reproducible RP-HPLC method for Simultaneous estimation of bulk and pharmaceutical formulations. Separation of Metformin and Linagliptin was successfully achieved on a THERMO, C18, 250cmx4.6mm, 5µm or equivalent in an isocratic mode utilizing KH_2PO_4 : Methanol (65:35) at a flow rate of 1.0mL/min and eluate was monitored at 226nm, with a retention time of 3.132 and 3.728 minutes for Metformin and Linagliptin respectively. The method was validated and found to be linear in the drug concentration range of 50µg/ml to 150 µg/ml for Metformin and 50µg/ml to 150 µg/ml for Linagliptin. The values of the correlation coefficient were found to be 0.999 for Metformin and 1 for Linagliptin respectively. The LOD and LOQ for Metformin were found to be 1.909 and 6.362 respectively. The LOD and LOQ for Linagliptin were found to be 0.0349 and 0.1163 respectively. This method was found to be good percentage recovery for Metformin and Linagliptin were found to be 100 and 100 respectively indicates that the proposed method is highly accurate. The method was extensively validated according to ICH guidelines for Linearity, Accuracy, Precision, Specificity and Robustness.

KEYWORDS: Metformin, Linagliptin, High performance liquid chromatography.

INTRODUCTION

Compound 1: Metformin is a oral tablet available as generic drugs and brand names are Glucophage, fortamet and glumetza. Metformin decreases hepatic glucose production, decreases intestinal absorption of glucose, and improves insulin sensitivity by increasing peripheral glucose uptake and utilization.

Compound 2: Brand name of drug is tradjenta and generic name is linagliptin. It is a DPP-4 inhibitor developed by Boehringer Ingelheim for the treatment of type II diabetes. Two pharmacological characteristics that sets linagliptin apart from other DPP-4 inhibitors is that it has a non-linear pharmacokinetic profile and is not primarily eliminated by the renal system.

MATERIALS AND METHOD

Chemicals and Reagent:

PREPARATION OF MOBILE PHASE

Transfer 1000 ml of HPLC water into 1000 ml of beaker and add 0.1M KH₂PO₄. Transfer the above prepared KH₂PO₄ buffer and Methanol is mixed in the proportion of (65:35). They are mixed and sonicated for 20 min.

PREPARATION OF METFORMIN AND LINAGLIPTIN STANDARD AND SAMPLE SOLUTION

PREPARATION OF STANDARD SOLUTION

Accurately weigh and transfer 500 mg Metformin and 20mg Linagliptin into 100 ml of volumetric flask and add 10 ml of methanol and sonicate 10 min (or) shake 5 min and make with methanol.

Transfer the above solution into 1 ml into 10 ml volumetric flask dilute to volume with water.

METHODOLOGY

PREPARATION OF SAMPLE STOCK SOLUTION

Commercially available 20 tablets were weighed, powdered and the powdered equivalent to the 870 mg of Metformin and Linagliptin of active ingredients were transferred into a 100 ml of volumetric flask and add 10 ml of Methanol and sonicate 20 min (or) shake 10 min and make up with methanol. Transfer above solution 1ml into 10 ml of the volumetric flask dilute the volume with Water. And the solution was filtered through 0.45µm filter before injecting into HPLC system.


RESULTS AND DISCUSSION

SYSTEM SUITABILITY: System suitability data of Metformin and Linagliptin

parameter	Metformin	Linagliptin	Acceptance criteria
Retention time	3.132	3.728	+/-10
Theoretical plates	4560	7688	>2500
Tailing factor	1.59	1.56	<2.00
% RSD	0.4	0.4	<2.00

Standard Results of Metformin

S. No.	Sample name	RT	Area	USP plate count	USP tailing
1.	Injection 1	3.728	4634008	7668	1.59
2.	Injection 2	3.729	4606703	7787	1.57
3.	Injection 3	3.726	4631981	7762	1.60
4.	Injection 4	3.723	4622848	7646	1.59
5.	Injection 5	3.724	4653812	7713	1.59

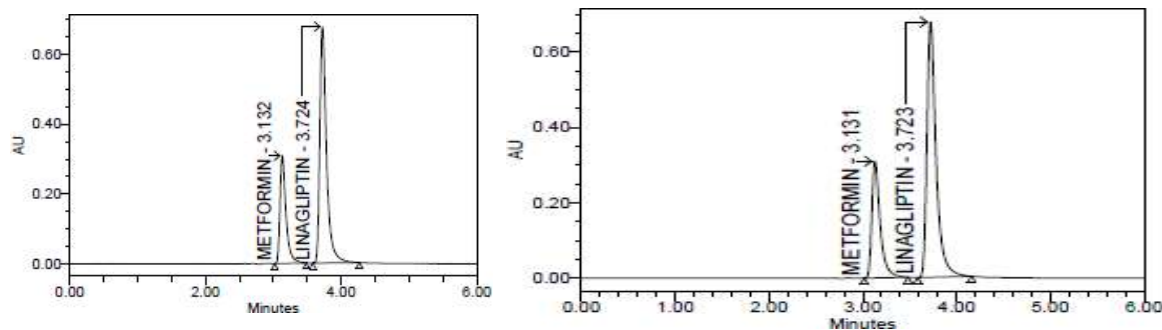

Typical Chromatogram of Standard-Typical Chromatogram of Sample-

RESULT

Results of system suitability study are summarized in the above table. Six consecutive injections of the standard solution showed uniform retention time, theoretical plate count, tailing factor and resolution for both the drugs which indicate a good system for analysis.

SPECIFICITY


S. No.	Sample name	Metformin Area	Rt	Linagliptin Area	Rt
1	Standard	1892041	1.132	4606966	3.728
2	Sample	1904192	3.131	4627816	3.723
3	Blank	-	-	-	-
4	Placebo	-	-	-	-

Results of forced degradation study for Metformin


Type of stress	Degradation products/Drug	Retention time	% Area	Peak purity	Result	% Assay	% Amount Degraded
Acidic Hydrolysis (mg/mL in 1N HCl) at 70°C for 2 days	-	3.130	1653385	0.999	Passed	86	14
Basic Hydrolysis (mg/mL in 1N NaOH) at 70°C for 2 days	-	3.130	1634097	0.999	Passed	85	15
Oxidative Hydrolysis (mg/mL in 3% v/v H ₂ O ₂) at 70 °C for 2 days	-	3.133	1643883	0.999	Passed	86	14
Photo Degradation (to UV light) for 14 days	-	3.131	1617526	0.999	Passed	84	16
Thermal Degradation at 70°C for 14 days	-	3.134	1608175	0.999	Passed	84	16

Results of forced degradation study for Linagliptin

Type of stress	Degradation products/ Drug (D)	Retention time	% Area	Peak purity	Result	% Assay	% Amount Degraded
Acidic Hydrolysis (mg/mL in 1N HCl) at 70°C for 2 days	-	3.726	3891067	0.999	passed	84	16
Basic Hydrolysis (mg/mL in 1N NaOH) at 70°C for 2 days	-	3.729	3911416	0.999	passed	84	16
Oxidative Hydrolysis (mg/mL in 3% v/v) at 70 °c for 2 days	-	3.731	3870137	0.999	passed	83	17
Photo Degradation (to UV light) for 14 days	-	3.730	3909913	0.999	passed	84	16
Thermal Degradation at 70°C for 14 days	-	3.733	3920769	0.999	passed	84	16

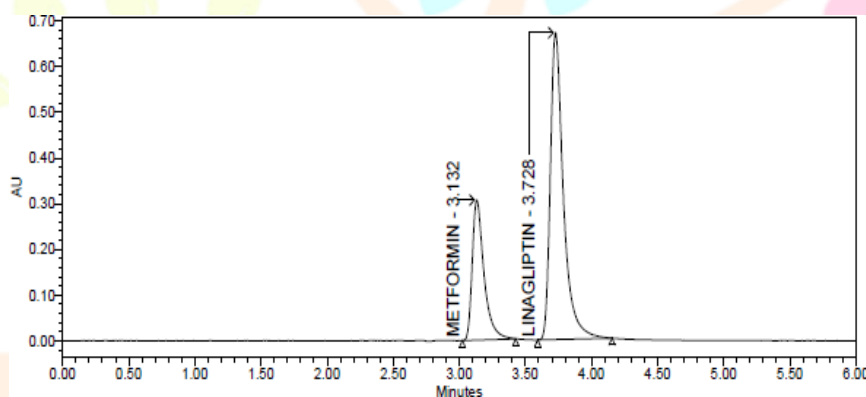

Chromatogram representing specificity of standard

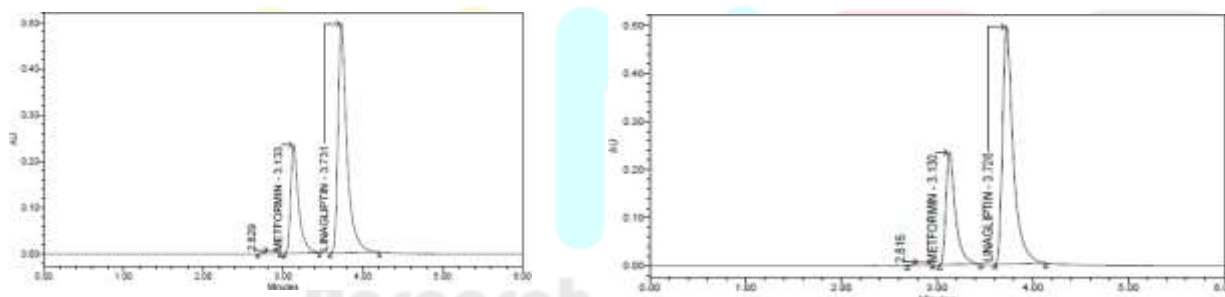

Chromatograms of Acid stress treated Metformin and Chromatograms of Base stress treated Metformin and Linagliptin mixture

RESULT

The forced degradation study showed the method was highly specific, the chromatographic peaks does not excipients have no effect on the analytical method. On the other hand, interfere with any other impurities. This proves that, blank peak did not overlap drug peak. blank peak did not overlap drug peak.

So the method is highly selectively.

ACCURACY :- Accuracy data for Metformin

S.NO	Accuracy Level	Injection	Sample area	RT
1	50%	1	953677	3.122
		2	953428	3.124
		3	953033	3.122
2	100%	1	1901769	3.131
		2	1901974	3.134
		3	1902392	3.136
3	150%	1	2868938	3.141
		2	2865114	3.152
		3	2860981	3.150

Accuracy (%recovery) results of Metformin


S.NO	Accuracy level	Sample name	Sample weight	µg/ml added	µg/ml found	% Recovery	% Mean
1	50%	1	435.00	247.500	245.91	99	99
		2	435.00	247.500	246.08	99	
		3	435.00	247.500	246.06	99	
2	100%	1	870.00	495.000	494.09	100	100
		2	870.00	495.000	492.55	100	
		3	870.00	495.000	493.48	100	
3	150%	1	1305.00	742.500	738.64	99	99
		2	1305.00	742.500	737.88	99	
		3	1305.00	742.500	736.69	99	

Accuracy data for Linagliptin

S.NO	Accuracy Level	Injection	Sample area	RT
1	50%	1	2325183	3.716
		2	2317701	3.716
		3	2317648	3.713
2	100%	1	4620300	3.721
		2	4626991	3.725
		3	4622070	3.726
3	150%	1	6948428	3.732
		2	6949946	3.744
		3	4940474	3.739

Accuracy (%recovery) results of Linagliptin

S.NO	Accuracy Level	Sample name	Sample weight	$\mu\text{g/ml}$ added	$\mu\text{g/ml}$ found	% Recovery	% Mean
1	50 %	1	435.00	10.000	9.90	99	99
		2	435.00	10.000	9.88	99	
		3	435.00	10.000	9.90	99	
2	100%	1	870.00	20.000	19.73	99	99
		2	870.00	20.000	19.74	99	
		3	870.00	20.000	19.75	99	
3	150%	1	1305.00	30.000	29.66	99	99

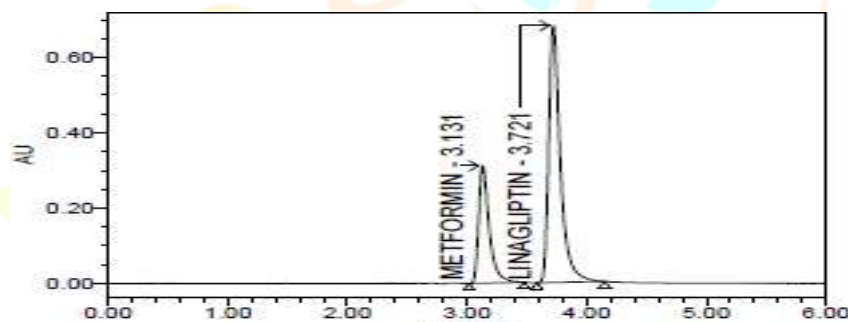
**Typical chromatogram for Accuracy 100 %****RESULT**

Results of accuracy study are presented in the abovetable. The measured value was obtained by recovery test. Spiked amount of both the drug were compared against the recovery amount. % Recovery was 100.00% for Metformin and 100.00% for Linagliptin. All the results indicate that the method is highly accurate.

PRECISION**Precision data for Metformin**

S. No	RT	Area	%Assay
injection1	3.131	1904192	99
injection2	3.132	19000711	99
injection3	3.137	1907020	99
injection4	3.134	1908231	99
injection5	3.127	1909733	99
injection6	3.131	1906386	99
Mean			99
Std. Dev.			0.17
% RSD			0.17

Precision data for Linagliptin

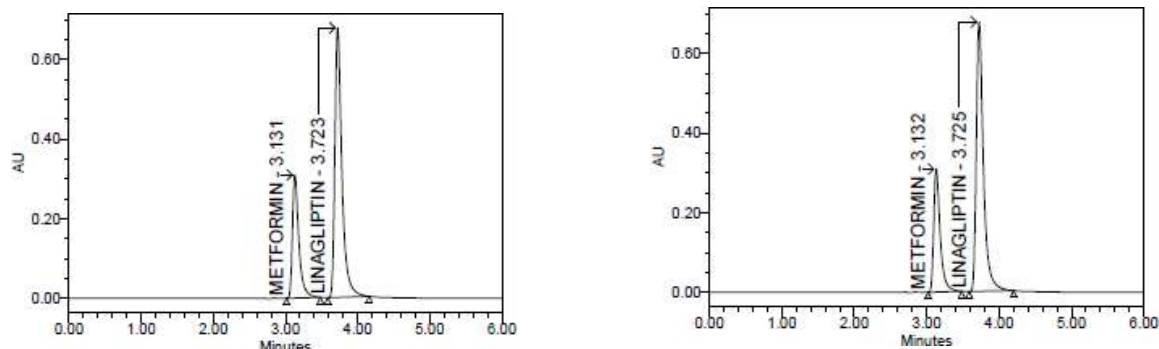

S. No	RT	Area	%Assay
injection 1	3.723	4627816	100
injection 2	3.725	4624364	100
injection 3	3.730	4628747	100
injection 4	3.725	4626814	100
injection 5	3.719	4626237	100
injection 6	3.721	4623058	100
Mean			100
Std. Dev.			0.05
%RSD			0.05

Chromatogram for precision injection

RESULTS



Results of variability were summarized run. Percentage relative standard deviation (%RSD) was found to be in the abovetable. % RSD of peak areas was calculated for various less than 2% which proves that method is precise.

LINEARITY

Linearity data for Metformin

S. No	Conc (µg/ml)	RT	Area
1.	50	3.125	953647
2.	75	3.137	1438138
3.	100	3.140	1902194
4.	125	3.145	2380153
5.	150	3.153	2867803
Correlation coefficient (r^2)			0.999

Linearity plot of Metformin

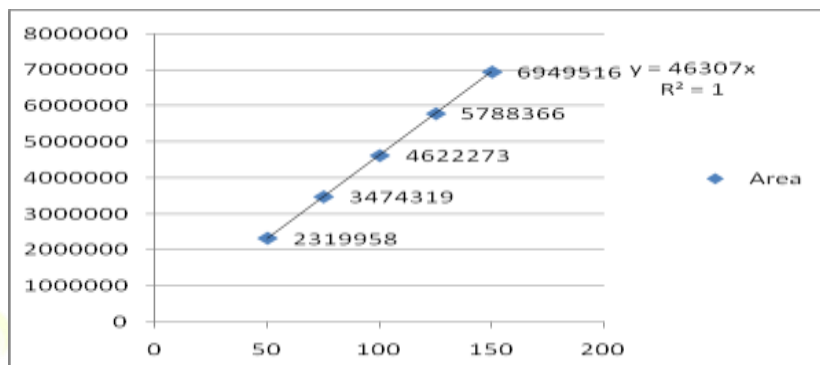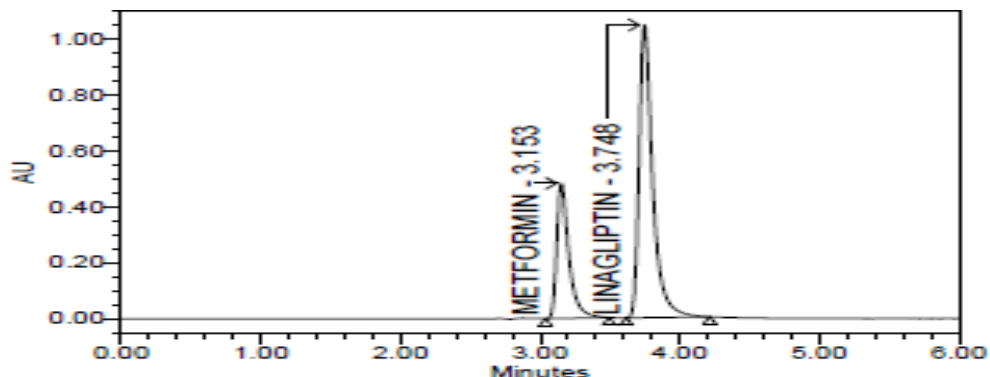

Chromatogram representing linearity

RESULT

A linear relationship between peak areas versus concentrations was observed for Glecaprevir and Linagliptin in the range of 50% to 150% of nominal concentration. Correlation coefficient was 0.999 for both Metformin and Linagliptin which prove that the method is linear in the range of 50% to 150%.


ROBUSTNESS:

Robustness data for Metformin


Parameter	RT	Theoretical plates	Asymmetry
Decreased flow rate(0.8ml/min)	3.945	7450	1.62
Increased flow rate(1.2ml/min)	2.621	6131	1.55
Decreased temperature(20 ⁰ c)	3.940	7434	1.61
Increased temperature(30 ⁰ c)	2.621	6131	1.52

Robustness data for Linagliptin


Parameter	RT	Theoretical plates	Asymmetry
Decreased flow rate (0.8ml/min)	4.678	8484	1.55
Increased flow rate (1.2ml/min)	3.118	7356	1.56
Decreased temperature(20 ⁰ c)	4.676	8409	1.55
Increased temperature(30 ⁰ c)	3.121	7304	1.55

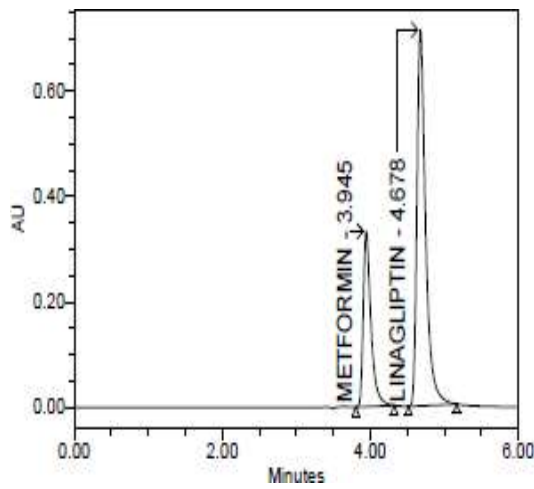


Chromatogram for decreased flow rate

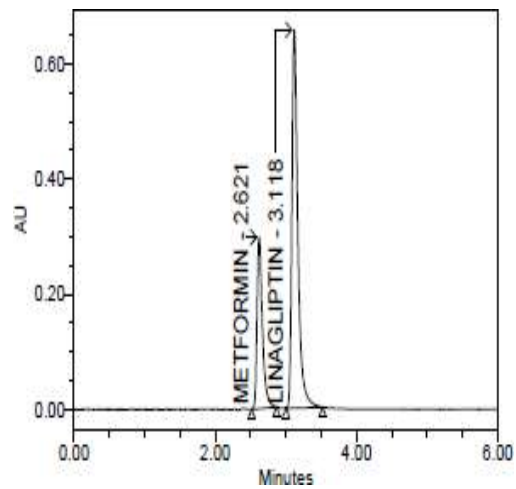

Chromatogram for increased flow rate

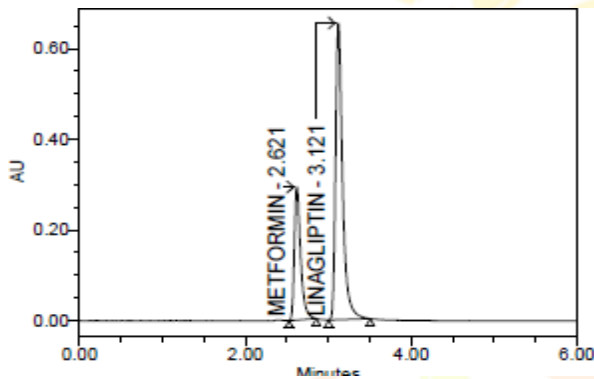

Chromatogram for decreased temperature

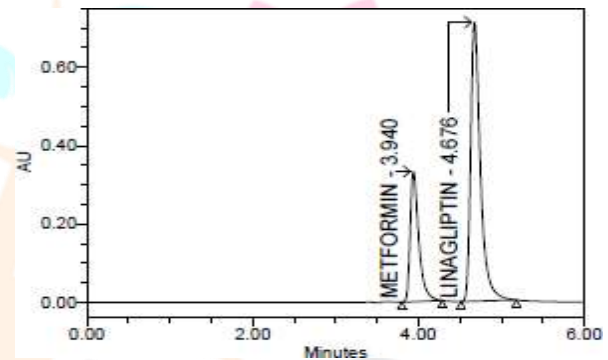

Chromatogram for increased temperature

RESULT

The results of Robustness of the present method had shown that changes made in the Flow and Temperature did not produce significant changes in analytical results which were presented in the above table. As the changes are not significant we can say that the method is Robust.

LIMIT OF DETECTION: Minimum concentration of standard component in which the peak of the standard gets merged with noise called the LOD $LOD = 3.3 * \sigma / S$

Where; σ = standard deviation S = slope LOD for Metformin = 1.909 LOD for Linagliptin = 0.0349

LOD data for Metformin and Linagliptin

S. No.	Sample name	RT	Area
1	Metformin	3.127	4887
2	Linagliptin	3.724	12240

Chromatogram for LOD:

LIMIT OF QUANTIFICATION

Minimum concentration of standard component in which the peak of the standard gets detected and quantification

$LOQ = 10 * \sigma / S$

Where; σ = standard deviation S = slope LOQ for Metformin = 6.362

LOQ for Linagliptin = 0.1163

LOQ data for Metformin and Linagliptin

S. No	Sample name	RT	Area
1	Metformin	3.134	3255
2	Linagliptin	3.716	9713

REFERENCES

1. H.H. Willard, L.L. Merritt, J.A. Dean, F.A. Settle. Instrumental Methods of Analysis, 7th edition, CBS publishers and Distributors, New Delhi, 1986; P.518-521, 580-610.
2. John Adamovics. Chromatographic Analysis of Pharmaceutical, Marcel Dekker Inc. New York, II Ed, P.74, 5-15.
3. Gurdeep Chatwal, Sahm K. Anand. Instrumental methods of Chemical Analysis, 5th edition, Himalaya publishing house, New Delhi, 2002; P.1.1-1.8,2.566-2.570.
4. D. A. Skoog. J. Holler, T.A. Nieman. Principle of Instrumental Analysis, 5th edition, Saunders College Publishing, 1998; P.778-787.
5. Skoog, Holler, Nieman. Principals of Instrumental Analysis, 5th Edition, Harcourt Publishers International Company, 2001; P.543-554.
6. William Kemp. Organic Spectroscopy, Palgrave, New York, 2005; P.7-10, 328-330.
7. P.D. Sethi. HPLC: Quantitative Analysis Pharmaceutical Formulations, CBS Publishers and distributors, New Delhi (India), 2001; P.3-137.
8. Michael E, Schartz IS, Krull. Analytical method development and Validation, 2004; P. 25-46.
9. R. Snyder, J. Kirkland, L. Glajch. Practical HPLC method development, II Ed, A Wiley International publication, 1997; P.235, 266-268, 351-353.653- 600.686-695.
10. Basic Education in Analytical Chemistry. Analytical Science, 2001: 17(1).
11. Method validation guidelines International Conference on harmonization; GENEVA, 1996.
12. Berry RI, Nash AR. Pharmaceutical Process Validation, Analytical method validation, Marcel Dekker Inc. New work, 1993; 57: 411-28.
13. Anthony C Moffat, M David Osselton, Brian Widdop. Clarke's Analysis of Drugs and Poisons, Pharmaceutical Press, London, 2004; PP 1109-1110,1601-1602.
14. Klaus Florey, Analysis Profile of Drugs Substances, Academic Press, New York, 2005; P.406-435.
15. P.N. Arora, P.K. Malhan. Biostatistics, Himalaya Publishers House, India, P.113, 139-140,154.
16. Doserge, Wilson and Gisvold's text book of organic medicinal and pharmaceutical chemistry, 8th edn, Lippincott Company, 1982; P.183-197.
17. British Pharmacopoeia, Government British Pharmacopoeial Commission, volume I and II, 2011, monograph 1188.

IJNRD
Research Through Innovation